

A

PATENT

Docket No. FWLPATUS012



Box Patent Application
Commissioner of Patents and Trademarks
Washington, D.C. 20231



NEW APPLICATION TRANSMITTAL
(STANDARD FORM)

Transmitted herewith for filing is the patent application of

Inventor(s): Binie V. Lipps, Frederick W. Lipps

For (title): Anti-LTNF for In Vitro Assay of Biological Toxins

1. **Type of Application**

This new application is for an Original Application.

2. **Papers Enclosed Which Are Required For Filing Date Under 37 CFR 1.53(b) (Regular) or 37 CFR 1.153 (Design) Application**

 17 Pages of specification

 3 Pages of claims

 1 Page of Abstract


 0 Sheets of drawings

CERTIFICATION OF EXPRESS MAILING DATE

I hereby certify that this correspondence is being deposited by me with the United States Postal Service on 27 April 1999 in an envelope as "Express Mail, Post Office to Addressee",

bearing Label Number EJ729250468US, addressed to the "Commissioner of Patents and Trademarks, Washington, D.C. 20231".

Date 4-27-99


John R. Casperson Reg. No. 28, 198

Send correspondence to:

John R. Casperson
PO Box 2174
Friendswood, Texas 77549

3. Additional papers enclosed

☒ PTO 1449s (3 sheets)
☒ Preliminary Amendment

4. Declaration or oath

☒ Enclosed
executed by
☒ inventors.

5. Language

☒ English

6. Small Entity Statement(s)

☒ Verified Statement that this is a filing by a small entity under 37 CFR 1.9 and 1.27 attached.


7. Fee Payment Being Made At This Time

☒ Enclosed

<input checked="" type="checkbox"/> basic filing fee	\$ 380.00
8 independent claims over 3 @ \$39.00	\$ 312.00
Total fees enclosed	\$ 692.00

8. Method of Payment of Fees

☒ A check in the amount of \$692.00 is attached.



John R. Casperson
Registration No. 28,198

Please send correspondence to:

John R. Casperson
PO Box 2174
Friendswood, Texas 77549

(281)-482-2961

PATENT

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

In re application of:
Binie V. Lipps
Frederick W. Lipps

Serial No.:

Filed:

For:

ANTI-LTNF FOR IN VITRO ASSAY
OF BIOLOGICAL TOXINS

§ ATTY DCKT NO: FWLPATUS012
§
§
§
§ Art Unit:
§
§ Examiner:
§
§
§
§

VERIFIED STATEMENT (DECLARATION) CLAIMING SMALL ENTITY
STATUS (37 CFR 1.9(f) AND 1.27(b))--INDEPENDENT INVENTOR

As a below named inventor, I hereby declare that I qualify as an independent inventor as defined in 37 CFR 1.9 (c) for purposes of paying reduced fees under Section 41(a) and (b) of Title 35, United States Code, to the Patent and Trademark Office with regard to the invention **entitled**:

ANTI-LTNF FOR IN VITRO ASSAY OF BIOLOGICAL TOXINS

described in

☒ the specification filed herewith

I have not assigned, granted, conveyed or licensed and am under no obligation under contract or law to assign, grant, convey or license, any rights in the invention to any person who could not be classified as an independent inventor under 37 CFR 1.9 (c) if that person had made the invention, or to any concern which would not qualify as a small business concern under 37 CFR 1.9(d) or a nonprofit organization under 37 CFR 1.9(e).

Each person, concern or organization to which I have assigned, granted, conveyed, or licensed or am under an obligation under contract or law to assign, grant, convey, or license any rights in the invention is listed below:

☒ no such person, concern, or organization

[] persons, concerns or organizations listed below

I acknowledge the duty to file, in this application or patent, notification of any change in status resulting in loss of entitlement to small entity status prior to paying, or at the time of paying, the earliest of the issue fee or any maintenance fee due after the date on which status as a small entity is no longer appropriate. (37 CFR 1.28(b)).

I hereby declare that all statements made herein of my own knowledge are true and that all statements made on information and belief are believed to be true; and further that these statements were made with the knowledge that willful false statements and the like so made are punishable by fine or imprisonment, or both, under Section 1001 of Title 18 of the United States Code, and that such willful false statements may jeopardize the validity of the application, any patent issuing thereon, or any patent to which this verified statement is directed.

Binie V. Lipps 4/22/99
Binie V. Lipps (date)

Frederick W. Lipps 4/22/99
Frederick W. Lipps (date)

PATENT

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

In re application of:	§	ATTY DCKT NO: FWLPATUS012
Binie V. Lipps	§	
Frederick W. Lipps	§	
Serial No.:	§	Art Unit:
Filed:	§	Examiner:
For:	§	
ANTI-LTNF FOR IN VITRO ASSAY	§	
OF BIOLOGICAL TOXINS	§	

PRELIMINARY AMENDMENT

Commissioner of Patents and Trademarks
Washington, D.C. 20231

Sir:

PRELIMINARY AMENDMENT

Prior to calculating the filing fee for the above identified patent application, kindly enter the following amendment.

IN THE SPECIFICATION

Kindly move the paragraph reading:

"The results of table 6 show that the detection level or the binding of C. atrox venom was 3.7 μ g in normal rabbit serum but dropped to 1.2 μ g in presence of homologous rabbit anti C. atrox serum. Thus, showing the difference of 2.5, which is the neutralizing index for this antiserum. The venom neutralized by the specific anti venom is not detected by anti-LTNF. The neutralizing index depends upon the potency of the anti serum. The neutralizing index for anti serum to V. russelli was 2.9 and for O. scutellatus 0.6"

from the top of page 17 to the bottom of page 16, so that the claims begin on a separate page.

IN THE CLAIMS

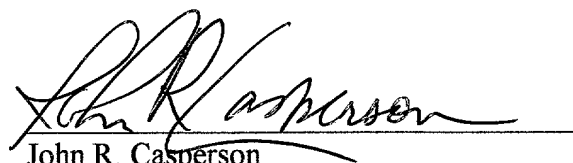
In claim 10, line 3 (page 18, line 30), delete "claims 1 and 2" and substitute therefor ---claim 2---.

In claim 16, line 1 (page 19, line 23), delete "Processes" and substitute therefor ---A process---

REMARKS

The amendment is made so that the application better complies with patent office requirements and more clearly does not include a multiple dependent claim.

Favorable consideration of the application at early date is respectfully requested.


John R. Casperson
Registration No. 28,198

Please send correspondence to:

John R. Casperson
P.O. Box 2174
Friendswood, Texas 77549

(281)-482-2961

APPLICATION FOR PATENT

INVENTORS: BINIE V. LIPPS,
FREDERICK W. LIPPS

5

TITLE: ANTI-LTNF FOR IN VITRO ASSAY
OF BIOLOGICAL TOXINS

BACKGROUND OF THE INVENTION

10

The invention relates to the introduction of a novel reagent for estimating biological toxins by in vitro test, as a replacement of currently used animal bioassay.

15 The direct testing of the toxins on animals poses ethical and moral questions. Therefore, there is great need for an alternative assay for biological toxins as a replacement for animal testing.

20 Biological toxins are grouped according to the source, such as: animal, plant, unicellular or one celled algae, and bacteria. Toxins are diverse and range from well-defined single macromolecules (tetanus, diphtheria, and botulinum toxins) to mixtures of complex molecules such as: snake or scorpion venoms or simple chemical entities (digitoxin, colchicine, tricyclic anti-depressants).

25

Mouse bioassay is the accepted and practiced method for assaying toxins, and is recommended by the Association of Official Analytical Chemists (AOAC). However, the mouse test for biological toxins is expensive and is disliked by animal activists and pressure is mounting worldwide to eliminate live animal bioassays. Germany, Switzerland and Australia have already banned the use of mammals for venom lethality tests.

30

The mouse bioassay suffers from other shortcomings as well. Toxins not lethal to mouse will go undetected by mouse bioassay. Also, mouse bioassay collectively detects the lethal effects of toxin, whether it is a single toxin or a mixture of several toxins.

5

Numerous investigators have developed immunological tests for assaying various toxins. Immunological tests are based on the principle of antigen-antibody reaction. Polyclonal antibodies are raised by immunizing animals: mouse, rabbit, rat, goat, etc. with a toxin.

10

Monoclonal antibodies are prepared in hybridoma cells. Spleen cells from an animal immunized with desired toxin are fused with corresponding species of myeloma cells. Both monoclonal and polyclonal antibodies produced are specific to the toxin (antigen used for immunization), and react with specific toxin.

15

The most frequently used in vitro immunological test is Enzyme Linked Immunosorbent Assay (ELISA). Biological toxins are assayed by several types of ELISA tests using specific antibodies to the desired toxin. ELISA tests can be carried out with numerous variations, although the most common format for detecting toxin from serum samples is the antigen-capturing method also known as double-sandwich method.

20

25

Li and Ownby reported the development of ELISA test for identification of venoms from snakes in the Agkistrodon genus. Several investigators have reported the development of ELISA for different toxins such as ricin and botulinum toxins; types A, B, and E in inoculated food samples. Morton and Tindall compared three tests; HPLC fluorescent method and two monoclonal antibody test kits for the determination of okadaic acid content of dinoflagellate cells and the results were not consistent. Since outbreaks of diuretic shellfish poisoning (DSP) may be caused by okadaic acid, methylokadaic acid, or a combination of these toxins, they concluded that both kinds of ELISA kits may underestimate total toxin effect in toxic shellfish.

30

35

All these methods suffer disadvantages, which have prevented their wide spread implementation particularly under the regulatory requirements. The mouse bioassay detects a wide range of known
5 and presumably unknown toxins. Therefore, it is unlikely that it will be abandoned completely in favor of existing bioassays, until an alternative is found which is similarly responsive.

What is needed is a reagent that can be used in a common protocol
10 for an in vitro assay test for biological toxins.

What is further needed is an in vitro test that can recognize all types of toxins, as a replacement of animal use.

15

OBJECTS OF THE INVENTION

An object of the present invention is to provide a novel reagent, specifically, Anti Lethal Toxin Neutralizing Factor (anti-LTNF), having
20 a property to bind to all types of toxins derived from animal, plant and bacteria and can be assayed by immunological in vitro test, specifically, ELISA.

Another object of the invention is to provide a method to use Anti-
25 LTNF having such a unique property to recognize all types of toxins in vitro test with sensitivity greater than in vivo mouse test.

Another object of the invention is to provide anti-LTNFs (made versus natural and synthetic LTNF) for use in a simple, cheap, ethically
30 accepted in vitro alternative as a replacement for the currently used animal bioassay.

Another object of the invention is to provide a method of identifying toxins by an ELISA sandwich utilizing anti-LTNF as primary antibody
35 and a set of secondary antibodies specific to the known toxins selected for the test.

SUMMARY OF THE INVENTION

In one embodiment of the invention, there is provided the novel composition of matter, anti-LTNF. Anti-LTNF is the antibody
5 produced biologically in response to natural Lethal Toxin neutralizing Factor (LTNF-n) or LTNF-s the synthetic active domain from LTNF-n. These materials are generically referred to herein as LTNFs. LTNF-n is a protein isolated from opossum serum having a molecular weight of 68 kDa and the amino acid sequence from the N-terminal of Leu Lys
10 Ala Met Asp Pro Thr Pro Pro Leu Trp Ile Lys Thr Glu. LTNF-s is a synthetic fragment of LTNF-n consisting of the portion of the N-terminal sequence of LTNF-n. A highly desirable active species of anti-LTNF is produced biologically in response to an LTNF-s having the first ten amino acids of the above sequence.

15 In another embodiment of the invention, an unknown toxin in fluid state is brought together with an anti-LTNF. The resulting immunological reaction indicates the presence of toxin in the fluid. The immunological reaction can determine the potency of the toxin.

20 Because anti-LTNF immunologically binds to all types of toxins, it can be used in vitro to test toxicity or to test for the presence of toxins as an alternative to animal testing. It can also be used for the identification of toxins by reacting with specific antitoxins under
25 investigation.

DETAILED DESCRIPTION OF THE INVENTION

Anti-LTNF is the antibody made versus LTNF-n and anti- LTNF-s is
30 made versus LTNF-s.

As used herein, LTNF-n refers to natural LTNF isolated from opossum serum having a molecular weight of about 68 kDa. The amino acid (AA) sequence from the N-terminal for the first fifteen amino acids of
35 LTNF-n is:

Leu Lys Ala Met Asp Pro Thr Pro Pro Leu Trp Ile Lys Thr Glu.

For additional details, see US Patent 5744449, the disclosure of which is incorporated by reference.

5 The synthetic LTNF comprising of 15 amino acids was identified as the active domain of the LTNF-n.

10 Synthetic LTNF or LTNF-s refers to a synthetic peptide having a molecular weight of less than 2,500 and at least a portion of the sequence set forth above. Preferably, LTNF-s has at least a 3-amino acid portion of the sequence, more preferably, at least a 5 amino acid portion of the sequence, and most preferably, at least a 7 amino acid portion of the sequence. Specific examples of LTNF-s include LTNF-15, LTNF-10 and LTNF-5. LTNF-15 is the 15 amino acid peptide having the first 15 amino acids of the above sequence and a molecular weight of 1,989 daltons. LTNF-10 is the 10 amino acid peptide having the first ten amino acids of above sequence with molecular weight of 1242 daltons. LTNF-5 is the 5 amino acid peptide having the first five amino acids of the above sequence with molecular weight of 647 daltons.

20

PRODUCTION OF ANTI-LTNF

25 Antibodies in general, are produced in animals by injecting the desired protein, which makes specific circulating antibodies to the protein. The serum of the immunized animal reacts immunologically to the protein, which was used for immunization. The presence of antibodies to the protein can be detected by various in vitro tests such as ELISA and immunoprecipitin test (IP).

30 LTNF-n and synthetic peptides consisting of fifteen, ten and five amino acids, designated as LTNF-15, LTNF-10 and LTNF-5 were used to produce antibodies by immunizing Balb/C mice. The sera of the immunized mice contained specific antibodies. In other words, mice immunized with LTNF-n produced antibodies specific to LTNF-n, anti
35 LTNF-n is formed. Likewise, mice immunized with synthetic peptides

LTNF-15, LTNF-10 and LTNF-5 produced anti-LTNF-15, anti-LTNF-10 and anti-LTNF-5 respectively.

- 5 (1) Polyvalent anti-LTNFs can be made in animals, mice, rabbits, goats etc. (2) monoclonal anti-LTNFs can be made in fused hybridoma cells secreting antibodies (3) or anti-LTNFs can be made in ascitic cavities of mice by injecting hybridoma antibody secreting cells.

10 There is a perception that small synthetic peptides do not generate antibodies on injection into animals. However, synthetic peptide can generate antibodies if it is tagged with a complete protein, before injecting to the animal. Landsteiner coined the term hapten for a low molecular weight, chemically defined compound which could induce antibody formation only when coupled to a larger carrier protein
15 molecule, before injecting. Thereby, the injected animal makes antibodies to both the hapten and the carrier protein.

Synthetic chemically defined LTNFs comprising of fifteen, ten and five amino acids can be considered as haptens and therefore, theoretically
20 should not induce antibodies if injected without carrier protein. However, surprisingly, synthetic LTNF consisting of fifteen amino acids was able to induce antibodies on injecting in mice. Encouraged by these results, mice were injected with synthetic LTNFs consisting of ten and five amino acids to induce antibodies, in a similar protocol.
25 Each synthetic LTNF consisting of fifteen, ten and five amino acids was able to induce specific antibodies as assayed by ELISA and IP.

SEPARATION OF ANTI-LTNF IgG

30 **Purification of IgG from Mouse anti-LTNFs:** Purification of immunoglobulin (IgG) can be accomplished by various methods and one of them is by high pressure liquid chromatography (HPLC). Anti sera were fractionated on HPLC from Tosoh Co. (Japan) using ionic exchange column from Polymer Laboratories (UK) and Trizma-HCl
35 gradient buffer pH 7.4. For a run, 20 mg of serum protein containing desired antibody was loaded for fractionation. The IgG fraction was

collected dialyzed and concentrated using apparatus from Spectrum Co. The concentrated fraction of IgG was refractionated on HPLC under identical conditions to obtain homogeneous preparation of IgG.

5

USES OF ANTI-LTNFs

Anti-LTNFs can be used to assay the toxicity of venom or venom toxins.

- 10 Anti-LTNFs can be used to assay the toxicity of toxins from sources other than venoms.

Anti-LTNFs can be used to detect the presence of multiple unrelated toxins.

15

Anti-LTNFs can be used to detect the presence of toxins in foods or biological fluids such as blood serum, saliva, urine, milk, tears etc.

20

Anti-LTNFs can be used to determine the potency of antivenoms and antitoxins.

25

The unique property of anti-LTNF to bind to all tested types of toxins immunologically in a fashion comparable to animal bioassay, is the theme of this patent. The protein LTNF-n and the synthetic peptides LTNF-15, LTNF-10 and LTNF-5 do not react immunologically with toxins by ELISA or IP tests. However, the anti-LTNFn, anti LTNF-15, anti LTNF-10 and anti LTNF-5 react immunologically with all tested toxins, as revealed by ELISA and IP tests. Moreover, the binding affinity of anti-LTNFn and the anti-synthetic LTNF-15, LTNF-10 and LTNF-5 roughly proportional to the

30

lethal dose of the toxin under testing.

35

Using the novel reagent anti-LTNF, ELISA tests can be performed in several different ways: Furthermore, anti-LTNFs can be used in crude form as anti serum, or it's purified IgG. However, to avoid background purified, it is ideal to use IgG for ELISA tests.

ELISA tests for assaying toxins can be performed by using different concentrations of toxins and reacting with constant amount of purified IgG from anti-LTNF.

5

ELISA Method in Antigen Capture Format: ELISA can be done in antigen capture format (also called double sandwich) where anti-LTNF IgG from one species (mouse) is used to capture the toxin and anti-LTNF from another species (rabbit) is used to recognize the toxin. For commercialization, two types of antibodies either polyclonal or

10

monoclonal, to LTNFs, natural, or synthetic, should be used. One for capturing the toxin and other for conjugating with alkaline phosphatase or peroxidase to recognize the toxin as in the antigen capture format.

15

Toxin detection from foods, sera or other body fluids must be done by antigen capture format. The second antibody can be conjugated directly to alkaline phosphatase or horse radish peroxidase to eliminate one step of washing the plate and to save one-hour incubation period. IgGs from all four anti-LTNFs to: LTNF-n, LTNF-15, LTNF-10 and LTNF-5 were tested in ELISA immunological tests and were found equally good.

20

25

30

35

EXAMPLES

Example I -- Production of polyclonal antibodies

5

Production of Polyclonal Antibodies to LTNFs (LTNF-15, LTNF-10 and LTNF-5): Adult Balb/C mice were used for immunization. The mice were used in compliance with the U.S. Public Health Service policy on humane care and use of animals. For the first injection each type of LTNF was mixed with equal amount of Freund's complete adjuvant and for the subsequent doses, the antigens were mixed with incomplete Freund's adjuvant. The mice were injected intramuscularly with 100 µg/mouse, four to six times, two weeks apart. At the end of the immunization period, mice were bled through ophthalmic vein and sera were separated.

10
15

Enzyme-Linked Immunosorbent Assay (ELISA) for Cross Reactivity of Anti-LTNFs:

ELISA tests were performed in 96 well microtiter plate. The plate was coated with one concentration of antigen. The wells of the plate were filled with 10 µg/ml in 0.05 M phosphate buffered saline (PBS) pH 7.3 and 100 µl/well, for each type of LTNFs such as natural LTNF (LTNF-n), synthetic peptide consisting of 15 amino acids (LTNF-15), consisting of 10 amino acids (LTNF-10) and consisting of 5 amino acids (LTNF-5). The plate was incubated for overnight at room temperature. After 18 to 24 hours the plate was washed three times (3X) with PBS and the plate was blocked with 0.25 ml/well of 3% Teleostean gelatin from cold water fish (Sigma) for 1/2 hour at RT.

20
25

Mouse anti-LTNFn, anti-LTNF-15, anti-LTNF-10 and anti-LTNF-5 diluted three fold were added to the coated plate and the plate was incubated at 37°C for 1 to 2 hours. The plate was washed 3X with PBS and horseradish peroxidase conjugated with mouse IgG made in goat (Sigma) was added and incubated for 1 hour at 37°C.

30
35

After which the plate was washed 3X with PBS and reacted with O Phenylendiamine Dihydrochloride (OPD) for color development. The test was read after 1/2 hour visually or preferably on ELISA plate reader. The results are presented in table 1.

5

Immunoprecipitation Test (IP):

IP was carried out as described by Ouchterlony. Agarose in 1-% concentration was dissolved by heating in normal saline and 7 ml of it was placed in 35-mm diameter petri dish. After the hardening of agarose, five wells were punched, one in the center and four peripheral wells, 1.0 cm apart from the central well. The central well was filled with one type of undiluted anti serum and the peripheral with four LTNFs having concentration 5 mg/ml. The wells were filled several times 50 µl each time, over the period of 72 hours. The results are presented in table 1.

15

Table 1: Immunological Cross Reactivity of Natural and Synthetic LTNFs by ELISA and by IP.

Antigen	ELISA titer				Immuno. Precipitation			
	Anti LTNF	Anti LT15	Anti LT10	Anti LT5	Anti LTNF	Anti LT15	Anti LT10	Anti LT5
Nat-LTNF	24300	600	2700	900	+	+	+	+
Syn-LTNF-15	8100	2700	2700	900	+	+	+	-
Syn-LTNF-10	8100	8100	8100	900	-	-	+	-
Syn-LTNF-5	8100	900	2700	600	-	-	+	-

20 The results of table 1 show that synthetic peptides consisting of as small as five amino acids are capable of producing antibodies in mice. The antibodies versus LTNF-15, LTNF-10 and LTNF-5 cross react with each other and with the natural LTNF-n by ELISA and by IP. Among

the three, anti-LTNF-10 proves to be better than anti-LTNF-15 and anti-LTNF-5.

- 5 It appears that the active domain of LTNF-n resides in ten amino acids. Therefore, LT-10 produced antibodies having high potency immunologically. LT-15 is may be too long having five amino acids not necessary for the active domain and LT-5 is too short of five amino acids, to generate antibodies having similar immunological potency.

10

Example II -- Immunological reaction of Anti-LTNF with Venom Toxin

- 15 **Immunological binding of anti-LTNF to venoms by ELISA:** The lethal dose was determined by injecting intraperitoneally 0.1 ml of venom in various concentrations in 20g ICR mice. ELISA test was carried in 96 well microplate. The wells of the microplate were coated with 0.1 ml of various concentrations of venom as antigen starting
- 20 from 100 μ g to 0.000564 (564 nanogram) diluted threefold in 0.05 M phosphate buffer saline pH 7.4 (PBS) and incubated for overnight at room temperature. After 18 to 24 hours the plate was washed three times (3X) with PBS and the plate was blocked with 0.25 ml/well of 3% Teleostean gelatin from cold water fish (Sigma) for 1/2 hour at RT.
- 25 The plate was washed 3X with PBS and 0.1 ml/well of 10 μ g/ml purified mouse anti-LTNF IgG was added. The plate was incubated at 37°C for 1 to 2 hours. And then, the plate was washed 3X with PBS and horseradish peroxidase conjugated with mouse IgG made in goat was added and incubated for 1 hour at 37°C. After which the plate
- 30 was washed 3X with PBS and reacted with O Phenylendiamine Dihydrochloride (OPD) for color development. The test was read after 1/2 hour visually or preferably on ELISA plate reader. The results are shown in table 2.

35

5

Table 2: Immunological binding of anti-LTNF-n to venoms.

Venom	Lethal dose μg	ELISA dose μg	Ratio L/E
<i>Crotalus atrox</i> (Rattlesnake)	300	3.7	81
<i>Naja.n. kaouthia</i> (Thailand cobra)	35	0.1	350
<i>Vipera russelli</i> (Common viper)	90	0.4	225
<i>Oxyuranus scutellatus</i> (Australian taipan)	3.5	0.02	175
<i>Androctonus australis</i> (Scorpion)	15	0.2	75
<i>Astrotia stokesii</i> (Sea snake)	4.0	0.03	133

10 The results of table 2 show that the toxicity of venoms was roughly proportional to the ELISA binding or ELISA titer. The mouse lethal dose for *C. atrox* venom was 300 μg and ELISA was 3.7 μg . Mouse lethal dose for sea snake *Astrotia stokesii* was 4 μg and its ELISA was 0.03 μg .

15

20

Example III -- Toxicity is proportional to binding affinity with toxins as well

5 Mouse Lethality of Toxins and Binding Affinity to Anti-LTNF IgG by ELISA: The lethal dose was determined by injecting intraperitoneally 0.1 ml of toxin at various concentrations in 20 g ICR mice. ELISA test was performed as described above using three fold concentrations of toxin and 10 µg/ml anti-LTNF IgG. Results are shown in table 3.

10

Table 3: Mouse lethality of toxins in micro grams and binding affinity to anti-LTNF IgG by ELISA

Toxin	Source	Lethal dose µg	ELISA dose µg	Ratio L/E
BoTx	<i>Clostridium botulinum</i>	1	0.04	25
Cholera	<i>Cholerae vibrio</i>	ND*	0.11	
Cobratoxin	<i>Naja n. kaouthia</i>	3	0.2	15
Crotoxin	<i>Crotalus d. terrificus</i>	5	1.1	4.5
Holothurin	<i>Actinopyga agrassizi</i>	200	1.8	111
Ricin	<i>Ricinus communis</i>	2	0.04	50
Taipoxin	<i>Oxyuranus scutellatus</i>	5	1.1	4.5
Toxic shock	<i>Staphylococcus aureus</i>	ND*	0.33	

* ND = not determined because these toxins are not lethal to adult mice.

15

The results of table 3 show that lethality of toxins is roughly proportional to the ELISA binding or ELISA titer. Mouse lethal dose for BoTx (botulinum toxin) was 1 µg and ELISA detection level was found to be 0.04 µg, whereas, for Holothurin mouse lethal dose was 200 µg and its detection level by anti-LTNF was 1.8 µg. Lower the ELISA binding means higher the detection level.

25

Example IV -- Multiple Toxins can be detected

Anti-LTNFs Detect the Presence of Multiple Unrelated Toxins:

5 Equal volumes of ricin at 10 $\mu\text{g}/\text{ml}$ was mixed with equal concentration of different toxins; BoTx, Cholera toxin and Staphylococcus aureus toxin. ELISA tests were carried using individual toxin and ricin plus other toxin as antigens. The toxins individually and the mixtures were diluted three fold. The binding was assayed by IgGs from anti-ricin and anti-LTNF. The results are shown
10 in table 4.

Table 4: Anti-LTNFs detect presence of several unrelated toxins

Toxin	Binding affinity, μg	
	Anti-Ricin	Anti-LTNF
Ricin	330	400
Ricin+BoTx	330	370
BoTx	0	37
Ricin+Cholera	330	110
Cholera	0	37
Ricin+Sta. aureus	330	12.3
Toxic shock	0	45

15 The results clearly show that anti-ricin detects only ricin by classical specific antigen-antibody reaction and non-of the other toxins mixed with it. However, anti-LTNF detects the presence of ricin as well as other toxins; BoTx, Cholera and Sta. aureus. Thus, anti-LTNF has the unique property of binding with all tested types of unrelated biological
20 toxins, individually or in mixture.

25

Example V -- Toxins in Food or animal body fluids can be detected

ELISA Test for Detection of Toxins from Foods by Anti-LTNF:

- 5 Human serum, human urine, chicken broth and dairy milk were spiked with 10 μ /ml of different toxins; BoTx, Ricin, Cholera and Sta. aureus. ELISA test was carried as described previously, using anti-LTNF IgG. The results are seen in table 5.

10 **Table 5: Anti-LTNF detects toxins from foods and body fluids**

Diluent	Detection Level in nanogram			
	BoTx	Ricin	Cholera	Sta. Aureus
PBS	333	333	111	333
Serum	11	11	11	11
Milk	220	110	111	111
Broth	330	110	111	111
Urine	330	220	12.3	37

BoTx: Botulium toxin from Clostridium botulinum.

Ricin: Plant toxin from Castor seeds.

Cholera: Toxin from Vibrio cholerae.

Sta. aureus: Toxin from Staphylococcus aureus.

- 15 The detection levels of BoTx, ricin, cholera and toxic shock was the lowest in human serum -- down to 11 ng (This corresponds to the highest sensitivity). Detection levels in milk were higher than observed in broth and urine (except for BoTx). The detection levels for
- 20 cholera and Sta. aureus were lower in urine than in milk and broth.

- Note: The concentrations of serum, milk, broth and urine diluted threefold in PBS for the obtained detection levels without added toxins were ELISA negative. In other words the positive reaction was not due
- 25 to the protein present in diluent but due to the presence of added toxin. Thus, anti-LTNFs are useful for in vitro assaying various types of food poisoning toxins from liquid and solid foods.

Example VI -- Antivenom and Antitoxin potency can be determined.

- Neutralization assay of Anti-venoms using anti-LTNF:** Currently, potency of antivenoms and anti toxins is measured by neutralization tests in mice. Neutralization test requires numerous mice because this test requires a determination of the LD₅₀ of venom alone and for the venom neutralized by the antiserum. At least four different concentrations of venom alone and a similar number or more for the venom to be neutralized by the antivenom must be used. In order to achieve statistical significance five to six mice are used per concentration. Thus, to determine the potency of a single antivenom fifty mice are required. (This is undesirable.)
- The protein concentration of each venom was adjusted to 200 µg/ml. Each venom was mixed with equal amount of normal rabbit serum and homologous respective rabbit anti serum. The mixtures were incubated at 37°C for one hour. The mixtures of venom with normal serum and specific anti serum were considered as antigens for ELISA. ELISA test was performed as described above. The results of detection levels of venom toxins in the mixtures are shown in table 6.

Table 6: Neutralization assay of venoms using anti-LTNF

Venom	Binding to anti-LTNF µg	Neutralization Index
<i>Crotalus atrox</i> + NS	3.7	2.5
<i>Crotalus atrox</i> + anti serum	1.2	
<i>Naja. n. kaouthia</i> + NS	1.2	1.16
<i>Naja. n. kaouthia</i> + anti serum	0.04	
<i>V. russelli</i> + NS	3.7	2.9
<i>V. russelli</i> + anti serum	0.8	
<i>O. scutellatus</i> + NS	0.8	0.6
<i>O. scutellatus</i> + anti serum	0.2	

NS= normal rabbit serum

- Homologous rabbit anti serum specific to each venom

The results of table 6 show that the detection level or the binding of C. atrox venom was 3.7 µg in normal rabbit serum but dropped to 1.2 µg in presence of homologous rabbit anti C. atrox serum. Thus, showing the difference of 2.5, which is the neutralizing index for this antiserum. The venom neutralized by the specific anti venom is not detected by anti-LTNF. The neutralizing index depends upon the potency of the anti serum. The neutralizing index for anti serum to V. russelli was 2.9 and for O. scutellatus 0.6

CLAIMS

We claim:

1. A composition of matter comprising Anti-LTNF, which is polyclonal or monoclonal antibody made versus:
natural LTNF-n having mol. wt. of 68 kilodalton and a partial amino acid sequence from the N-terminal of fifteen amino acids
Leu Lys Ala Met Asp Pro Thr Pro Pro Leu Trp Ile Lys Thr Glu.
2. The composition of matter comprising Anti-LTNF including polyclonal or monoclonal antibodies made versus
any active portion of LTNF-n sequence and specifically those antibodies to LTNF-15, LTNF-10 and LTNF-5, comprising of 15, 10 and 5 amino acids, respectively, from the N-terminal of LTNF-n.
3. A composition of matter as in claim 2, wherein the polyclonal anti-LTNF is further characterized as immunoglobulin (IgG) from immunized animals and monoclonal anti-LTNF in claim 2 is further characterized as immunoglobulin (IgG) from fused hybridoma cells from immunized mouse or mouse myeloma cells, or other appropriate cells.
4. The in vitro toxin assay process, based on the use of Anti-LTNF (made versus LTNF-n, LTNF-15, LTNF-10 or LTNF-5) as a reagent for the detection of biological toxins, becomes an ethical replacement for

currently used animal bioassay, typically mice, or other rodents and other animals.

5 **5.** The anti-LTNFs made versus natural LTNF and versus synthetic peptides consisting of at least five amino acids detect biological toxins derived from animal, plant and bacteria by ELISA.

10 **6.** The anti-LTNFs made versus LTNF-n and the synthetic peptides comprising of 15, 10 and 5 amino acids recognize toxins derived from animal, plant and bacteria by ELISA assay.

15 **7.** The anti-LTNFs provide essential reagents for the in vitro assay of the wholesomeness of toxins existing in singular form, or in mixture, in a manner comparable to animal bioassay.

8. The anti-LTNFs detect and assay the toxins from foods, blood sera and other body fluids saliva, milk, urine etc. by ELISA test in antigen capture format, or any similar test.

20 **9.** The neutralizing potency of an anti-toxin is the neutralizing index given by the toxin assay minus an anti-toxin mixture assay; wherein, the toxin assay is determined by ELISA test of the toxin plus normal serum; and the anti-toxin mixture assay is determined by ELISA test of a mixture of toxin plus anti-toxin mixture, such mixture
25 containing a reduced amount of free toxin due to neutralization by the anti-toxin.

10. The neutralizing potency of anti-toxins including anti-venoms can be assayed by in vitro test using anti-LTNF compositions as in
30 claims 1 and 2 as reagent, and thus saving thousands of mice as well as time and money.

11. A composition of matter comprising an antibody made versus a peptide containing at least five amino acids from the N-terminal
35 sequence

Leu Lys Ala Met Asp Pro Thr Pro Pro Leu Trp Ile Lys Thr Glu.

12. A composition of matter as in claim 11, which is in the form of
an immunoglobulin selected from the group consisting of an
immunized animal serum, a hybridoma cell culture and a mouse
5 ascitic fluid.

13. A composition of matter as in claim 12, which reacts
immunologically with a toxin selected from the group consisting of an
animal toxin, a plant toxin and bacterial toxin.

10

14. A process comprising contacting, in vitro, a biological toxin with
an antibody made versus a sequence of at least five amino acids from
the N-terminal of the sequence

Leu Lys Ala Met Asp Pro Thr Pro Pro Leu Trp Ile Lys Thr Glu

15 under conditions to cause the biological toxin to react
immunologically with said antibody.

15. A process as in claim 14, wherein the novel antibody is made
against LTNF having a non-immunological binding with toxins such
20 that its antibody has the property of being able to react
immunologically in vitro with a wide range of biological toxins.

16. Processes as in claim 15 which is carried out according to an
ELISA double-sandwich method protocol.

25

17. A process as in claim 15, wherein the biological toxin reacts
immunologically with said novel antibody to produce a first reaction
product, said process further comprising contacting the first reaction
product with a set of antibodies made for specific known toxins for the
30 purpose of producing a modified ELISA test capable of identifying
toxins with the help of the novel anti-LTNF.

35

ABSTRACT OF THE DISCLOSURE

Lethal Toxin Neutralizing Factor has been isolated in purity from
5 opossum serum by high pressure liquid chromatography (HPLC)
fractionation. The amino acid sequence from the N-terminal for the
first fifteen amino acids of LTNF-n is: Leu Lys Ala Met Asp Pro Thr Pro
Pro Leu Trp Ile Lys Thr Glu. Antibodies to LTNF-n and synthetic
10 peptides consisting of fifteen, ten and five amino acids from the N-
terminal of the above sequence, designated as LTNF-15, LTNF-10 and
LTNF-5 were produced by immunizing Balb/C mice to produce Anti-
LTNF-n, Anti-LTNF-15, Anti-LTNF-10 and Anti-LTNF-5. The anti
LTNF-n, anti-LTNF-15, anti-LTNF-10 and anti-LTNF-5 react
15 immunologically with all types of toxins derived from animal, plant
and bacteria and can be assayed by immunological in vitro test such
as ELISA tests. Anti-LTNFs react roughly proportional to lethal dose
of biological toxins under in vitro immunological ELISA test similar to
the mouse bioassay test. This property of anti-LTNFs can be utilized
20 to replace the use of animals for bioassay of toxins from animal, plant
and bacteria. Anti-LTNFs can be polyclonal raised in animals or
monoclonal made in hybridomas. Anti-LTNFs can be used in crude
form for immunological in vitro testing. However, purified IgG from
the anti LTNFs is most desirable for consistent results from the in
vitro tests. In general, the in vitro use of Anti-LTNFs is a replacement
25 for animal use as is currently required for the assay of biological
toxins.

30

35

COMBINED DECLARATION AND POWER OF ATTORNEY
(JOINT INVENTORS)

As a below named inventor, I hereby declare that:

TYPE OF DECLARATION

This declaration is of the following type:

☒ original

INVENTOR IDENTIFICATION

My residence, post office address and citizenship are as stated below next to my name. I believe I am an original, first and joint inventor of the subject matter which is claimed and for which a patent is sought on the invention entitled:

TITLE OF INVENTION

ANTI-LTNF FOR IN VITRO ASSAY OF BIOLOGICAL TOXINS

SPECIFICATION IDENTIFICATION

the specification for which is attached hereto.

ACKNOWLEDGMENT OF REVIEW OF PAPERS AND DUTY OF CANDOR

I hereby state that I have reviewed and understand the contents of the above identified specification, including the claims, as amended by any amendment referred to above.

I acknowledge the duty to disclose information

- which is material to patentability as defined in 37, Code of Federal Regulations, § 1.56.

POWER OF ATTORNEY

I hereby appoint the following attorney(s) and/or agent(s) to prosecute this application and transact all business in the Patent and Trademark Office connected therewith:

John R. Casperson, Reg. No. 28,198.

SEND CORRESPONDENCE TO

John R. Casperson
PO Box 2174
Friendswood, Texas 77549

DIRECT TELEPHONE CALLS TO:

John R. Casperson - (281) 482-2961

DECLARATION

I hereby declare that all statements made herein of my own knowledge are true and that all statements made on information and belief are believed to be true; and further that these statements were made with the knowledge that willful false statements and the like so made are punishable by fine or imprisonment, or both, under Section 1001 of Title 18 of the United States Code, and that such willful false statements may jeopardize the validity of the application or any patent issued thereon.

SIGNATURE(S)

Full name of first inventor:

BINIE

V.

LIPPS

Given Name

Middle Initial or Name

Last Name

Inventor's signature

Binie V. Lipps

Date: April 22, 1999 Country of Citizenship: U.S.A.

Residence: 4509 Mimosa Drive, Bellaire, Texas 77401

Post Office Address: 4509 Mimosa Drive, Bellaire, Texas 77401

Full name of second inventor:

FREDERICK

W.

LIPPS

Given Name

Middle Initial or Name

Last Name

Inventor's signature

Frederick W. Lipps

Date: April 22, 1999 Country of Citizenship: U.S.A.

Residence: 4509 Mimosa Drive, Bellaire, Texas 77401

Post Office Address: 4509 Mimosa Drive, Bellaire, Texas 77401